




Roger Table Mic II is a microphone dedicated for working adults who participate in various meetings. It selects the person who's talking and switches automatically. Multiple Roger table microphones can be connected to create a network, making it ideal for large meeting configurations. It can also transmit the sound of multimedia e.g. computer.

Description

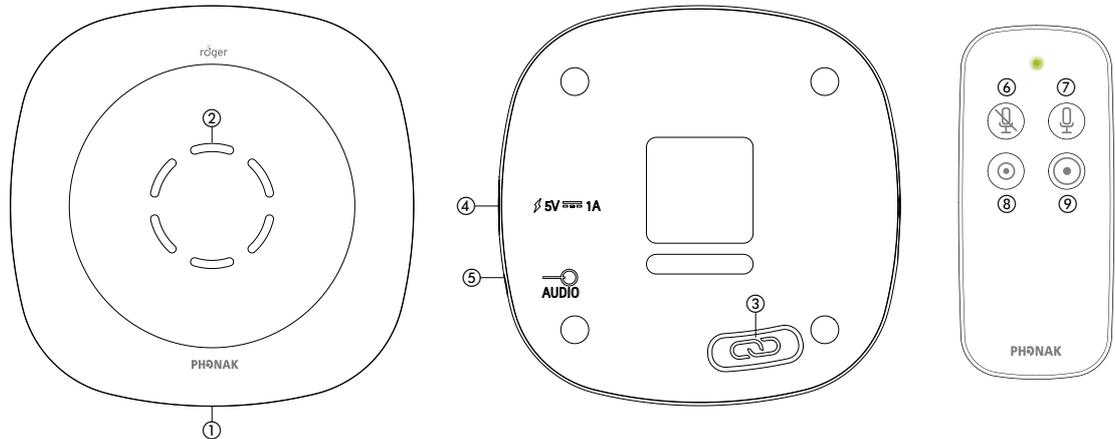

- 1) On/off button
- 2) Microphones and indicator lights
- 3) Connect button
- 4) Micro-USB socket
- 5) Audio input
- 6) Mute on
- 7) Mute off
- 8) Short pick-up range
- 9) Long pick-up range

Roger Table Mic II features

- MultiBeam Technology for improved speech in noise
- MultiTalker Network for easy coverage of large meetings
- Remote control for muting and adjusting pick-up range
- Audio input for connection to multimedia devices
- 16 hours operating time

General information

Dimensions:	69 x 69 x 14 mm / 2.7" x 2.7" x 0.6"
Weight:	103 g / 3.63 oz
Color:	Champagne
Operating conditions:	0°C to 45°C (+32°F to +113°F). Relative humidity of <90% (non-condensing)
Transport and storage conditions:	-20°C to 60°C (-4°F to +140°F). Relative humidity of <90% for a long period of time

Power supply

Type:	DC, stabilized
Primary voltage:	100 - 240 VAC
Secondary voltage:	5 VDC
Max. current:	1.0 A
Plug:	micro-USB
Important notice:	Only use the original Phonak power supply or an USB certified power supply 5 VDC, > 500 mA

Battery

Type:	Lithium Polymer
Capacity:	1000 mAh
Voltage:	3.7 V
Operating time:	16 hours

Audio information

Audio input:	Analog
Audio bandwidth:	100 Hz – 7.3 kHz

Microphone characteristics

Number of Microphones:	3
Microphone mode:	Table mode

Radio information

Streaming technology:	Roger technology
Frequency range:	2.4 GHz ISM band including adaptive automatic frequency hopping
RF power:	70 mW max.
Transmission range:	20 meters / 66 feet (typical) 40 meters / 130 feet (ideal conditions)
Max. paired devices:	Unlimited

Receiver compatibility

Compatible Roger receivers:	All Roger receivers
Max. connected receivers:	Unlimited

Microphone compatibility (MultiTalker Network)

Compatible microphones:	Roger Table Mic II Roger Select™ Roger Pen™ Roger EasyPen Roger Clip-On Mic
Max. connected microphones:	Up to 10

Standards

Europe	EN 300 328 EN 301 489 EN 60950-1 IEC/EN 62479
Canada	RSS-247 RSS-102
Japan	ARIB-T66
USA	CFR 47, part 15.231, part 15.247 FCC OET Bulletin No.65
