

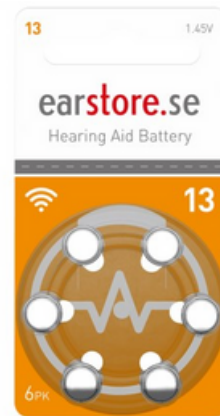
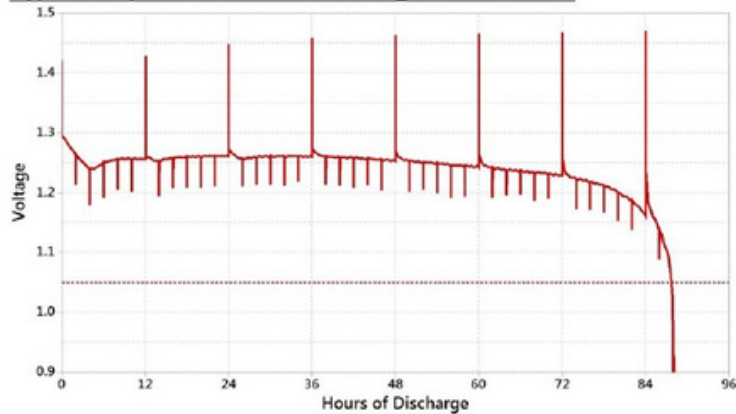


Typical IEC/ANSI Standard Discharge Performance

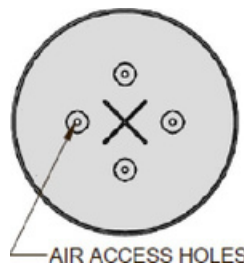

Article number	ES-13
Name	Earstore size 13 Orange Hearing Aid Battery
Cell Chemistry	Zinc Air
Recommended Use	Hearing Aid Instruments
Nominal Voltage	1.45 volts
Designation	IEC : PR48 ANSI : 7000ZD Common : 13
Weight	0.80 grams (0.028 oz)
Volume	0.26 cc (0.016 cu. in.)
Operating Temperature Range	-10° C to 50° C (14° F to 122° F)
Impedance	7.5 ohm @ AC 1000 Hz
Hearing Aid	144 hours – 2 mA background, 6 mA 100 ms pulse, once every 2 h, 12 h/d to 1.05V at 21°C(70°F) and 50% R.H.
IEC/ANSI Standard	86 hours – 3 mA background, 12 mA 100 ms pulse, once every 2 h, 12 h/d to 1.05V at 21°C(70°F) and 50% R.H.
IEC/ANSI Wireless	65 hours – 3 mA background 45 min, 5 mA 15 min pulse, every h, 12 h/d to 1.1V at 21°C(70°F) and 50% R.H.
Storage Guarantee	Three years with tab applied when stored at 21°C (70°F).
Activation and use	Remove and discard the plastic tab covering air holes. Allow the battery to stand for one minute to insure proper activation.
Environmental	These cells have no added mercury.
Compliant to the Following Standards	IEC 60086-1, IEC 60086-2 ANSI C18.1M Parts 1 & 2

Zinc Air Cell

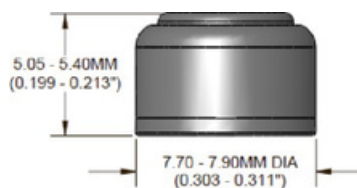

MDR 2017/745
Class 1

This data is subject to change. Performance information is typical. Contents herein do not constitute a warranty.

Accessory Line Int AB | Grafikgatan 8 | 754 50 Uppsala | +46(0)854553630

www.earstore.se